




CMOS+LED

FLEXIBLE VIDEO BRONCHOSCOPE



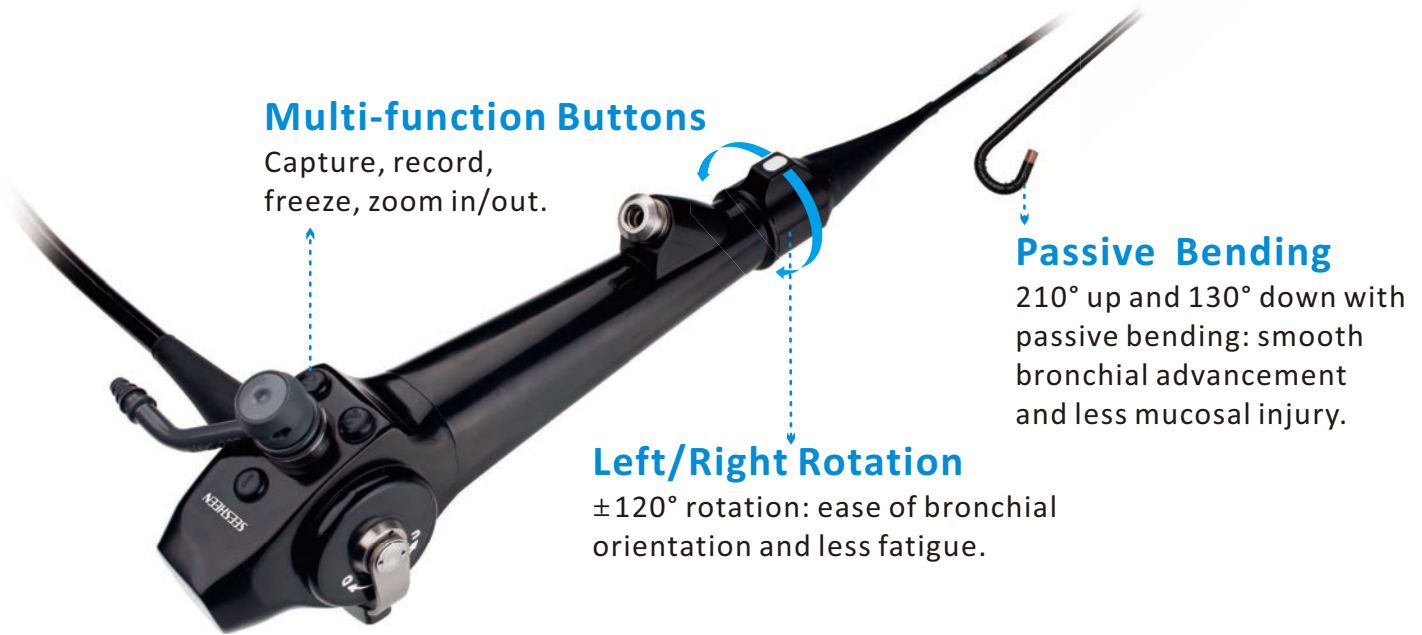

Multi-function Buttons

Capture, record, freeze, zoom in/out.

Passive Bending

210° up and 130° down with passive bending: smooth bronchial advancement and less mucosal injury.

Left/Right Rotation

±120° rotation: ease of bronchial orientation and less fatigue.

Portable

Compatible with portable monitor CV-1080/HV-3101.



CMOS + LED Technology

CMOS image sensor delivers a high quality image. No fiber: durable and easy maintenance.



Application

Pneumology: bronchoscopy, peripheral observation and interventional therapy.



Model	Insertion Tube(mm)	Channel (mm)	Working Length(mm)	Bending (up/down)	Field of View	Depth of Field(mm)	Rotation	Image Processor
BR-1231	3.1	1.2	600	210°/130°	120°	3-50	L120°/R120°	HV-3101/ CV-1080
BR-1242	4.2	2.0	600	210°/130°	120°	3-50	L120°/R120°	HV-3101/ CV-1080
BR-1249	4.8	2.0	600	210°/130°	90°	3-50	L120°/R120°	HV-3101/ CV-1080
BR-1259	5.9	2.8	600	180°/130°	90°	3-50	L120°/R120°	HV-3101/ CV-1080